





Changed
Cartridge System

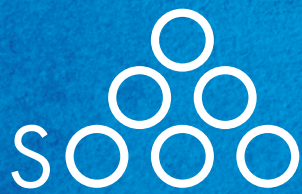

EFFICIENT IS **SILGINAT**[®]

ANATOMIC IMPRESSION

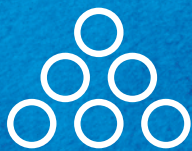

KETTENBACHDENTAL

Simply intelligent

SILGINAT® – SIMPLE, FAST, EFFICIENT!

Anatomical impression

Silginat® is an alginate replacement material based on A-silicone especially for **alginate applications**.









INDICATIONS

- Opposing dentitions
- Fabricating temporary crowns and bridges
- Anatomical models
- Fabricating simple removable prosthetic restorations
- Producing models for the construction of trays for whitening solutions, bites and surgical guides
- Orthodontic work
- Case study models
- Fabricating removable retainers and splints



Technical data

						
	Working time at 23 °C / 74 °F ≤	Intraoral setting time at 35°C / 95 °F ≥	Total setting time* ≥	Elastic recovery test (approx.) %	Linear dimensional change (approx.) %	Hardness (ap- prox.) Shore
 Silginat® / Silginat® Strawberry 1:1 cartridge	1:15 min.	1:15 min.	2:30 min.	99.5	- 0.20	A 42
 Silginat® / Silginat® Strawberry 5:1 cartridge	1:30 min.	1:30 min.	3:00 min.	99.5	- 0.20	A 45

*Total setting time (removal time from mouth) from the start of mixing

AN A-SILICONE AS YOU KNOW IT



- Storage stability and dimensional stability, high precision

Silginat vs. alginate
approx. 24 h after impression-taking



EASY TO HANDLE



- Clean and fast work thanks to the 5:1 jumbo or 1:1 50 mL cartridge with reproducible results



COST-EFFECTIVE IMPRESSION TAKING

- Pourable multiple times, only one anatomical impression per patient required



FOR EVEN GREATER EFFICIENCY

Saves time and many work steps during your work!



SILGINAT® – THE RANGE

Silginat®

Refill pack	2 x 380 mL cartridges	REF 14713
Normal pack	6 x 50 mL cartridges, 6 mixing tips	REF 1384611
Bulk pack	24 x 50 mL cartridges	REF 1384711

Silginat® Strawberry

Refill pack	2 x 380 mL cartridges	REF 14715
Normal pack	6 x 50 mL cartridges, 6 mixing tips	REF 1382611
Bulk pack	24 x 50 mL cartridges	REF 1382711

Accessories

Mixing tips blue, Ø 6.0 mm	50 tips	REF 17244
Mixing tips blue, Ø 6.0 mm	100 tips	REF 17246
Dynamic mixers white	45 tips	REF 17902

Changed
Cartridge System



NOTES

Kettenbach LP
16052 Beach Blvd. Ste. 221
Huntington Beach, CA 92647, USA
Phone: 877-532-2123
Fax: 888-892-9820
Website: www.kettenbach-dental.us
Web shop: www.kettenbachusa.com
Custservice@kettenbachusa.com

Not for use, sale and
distribution in Germany.

